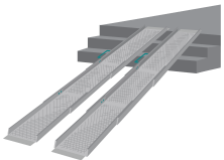
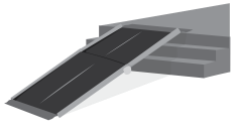

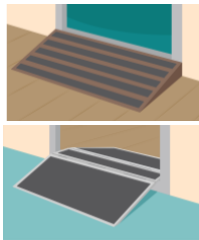



## Selection table: Portable ramps

| <b>Step 2: Select</b><br>Always select a portable ramp with the person.<br>If selecting a ramp for a public building or public transport, consider all the mobility assistive products available. |   | <b>Twin track portable ramp</b><br> | <b>Wide platform portable ramp</b><br> | <b>Permanent ramp</b><br> | <b>Other ramps</b><br> | <b>Other option for access</b><br> |
|---|---|---|---|--|---|---|
| <b>Environment</b>  | Ramp used in home                               | ✓   | ✓   | ✓  | ✓   | ✓   |
|   | Ramp used for public access                     | ✗   | ✓   | ✓  | ✓   | ✓   |
| <b>Doorways</b>   | Step on one side                                | ✓   | ✓   | ✓  | Consider wedge ramp   | Consider grab bars  |
|   | Step on both sides                              | ✗   | ✗   | ✓  | Consider two wedge ramps or hinge ramp  |   |
| <b>Surfaces</b>   | Surfaces uneven and/or not solid in all weather | ✗   | ✗   | ✓  | ✗   | Discuss ability to change surface   |
| <b>Measurements</b>   | Height more than 45cm                           | ✗   | ✗   | ✓  | ✗   | ✓   |
|   | Length more than 350cm                          | ✗   | ✗   | ✓  | ✗   | ✓   |
| <b>Mobility assistive products</b>  | Walking aid                                     | ✗   | ✓   | ✓  | ✓   | Consider grab bars  |
|   | Three-wheel wheelchair                          | ✗   | ✓   | ✓  | ✓   | ✓   |
|   | Four-wheel wheelchair                           | ✓   | ✓   | ✓  | ✓   | ✓   |
| <b>Lifestyle</b>  | Person will be assisted on ramp                 | ✗   | ✓   | ✓  | ✓   | ✓   |
|   | Person <b>cannot</b> lift and carry ramp        | ✗   | ✗   | ✓  | ✓   | ✓   |
|   | Person <b>cannot</b> store ramp                 | ✗   | ✗   | ✓  | ✓   | ✓   |