



1. Information about the person			
Family name ✎		Given names ✎	
Age: 0-5 ☐ 6-18 ☐ 19-39 ☐ 40-54 ☐ 55+ ☐		Gender: Male ☐ Female ☐ Other ☐	
Address ✎		Telephone ✎	

2. Assessment

Environment

Where will the ramp be used?	Entrance to: Home ☐ Personal vehicle ☐ Public building ☐ Public transport ☐
	Inside: Home ☐ Public building ☐
	Other: ✎

Will the ramp be positioned in a doorway?	Yes ☐ No ☐	Yes → Check width of mobility assistive product and ramp is less than width of doorway. Discuss other options for access if needed. No → Continue.
Is the supporting surface solid? <i>Consider surface in all weather conditions.</i>	Yes ☐ No ☐	Yes to both → Portable ramp may be suitable. No to either → Discuss if different surface is possible. If not possible, consider permanent ramp.
Are the supporting surfaces even?	Yes ☐ No ☐	

Measure

Height: <i>Measure the total height of step/s using a metal tape measure. Measurement in centimetres (cm).</i>	Height (in cm) = ✎	Is total height of the steps 45 cm or less? Yes ☐ No ☐	Yes → Continue. No → Consider permanent ramp.
Length: <i>Work out lengths for different ratios. Then circle lengths of 350 cm or less.</i>	1:8 Height x8 = ✎	1:10 Height x10 = ✎	1:12 Height x12 = ✎

Mobility assistive products

What mobility assistive product would you like to use on the ramp? <i>Select all that apply.</i>	Walking aid ☐	→ A wide platform portable ramp may be suitable.
	Three-wheel wheelchair ☐	
	Four-wheel wheelchair ☐	→ A wide platform or twin track portable ramp may be suitable.

Activity and ability *Ask the person the following questions.*

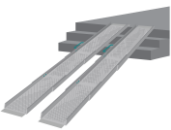
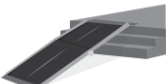
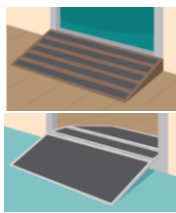
Do you need assistance on the portable ramp?	Yes ☐ No ☐	Yes and assistance is available → Continue. Consider a wide platform portable ramp. Yes and assistance is not available → ☹ Stop assessment and discuss other options for access. No → Continue. Consider a wide platform or twin track portable ramp.
Are you and/or an assistant able to carry and position a portable ramp?	Yes ☐ No ☐	Yes → A portable ramp may be suitable. No → Discuss other options for access.
Are you able to store the selected portable ramp?	Yes ☐ No ☐	Yes → A portable ramp may be suitable. No → Discuss other options for access.

Weight on ramp <i>Check common weights and weight limits at end of form.</i>			
Measure: <i>Measurements in kilograms (kg)</i>	(A) Weight of person Weight of assistant	A = kg	Is the total weight on ramp within weight limit of ramp? Yes ☐ → Continue. No ☐ → Consider permanent ramp.
	(B) Weight of mobility assistive product	B = kg	
	(A+B) Total weight on ramp	A+B = kg	

4. Plan			
Provide	Portable ramp: Twin track ☐ Wide platform ☐	Length:	Weight limit:
Teach	Lift and carry ☐ Position ☐	Discuss	Adjusting surface ☐ Other options ☐
Refer	Permanent ramp ☐ Wedge or hinged ramp ☐ Modify threshold ☐ Grab bars ☐ Other:	Follow up date	

Common weight of mobility assistive products	
Walking aid	Up to 10 kg
Manual wheelchair	8-25 kg
Electric wheelchair	60-130 kg

Common weight limits of portable ramps	
Wide platform portable ramp	300 kg
Twin track portable ramp	250-350 kg

Selection table: Portable ramps						
Step 2: Select Always select a portable ramp with the person. If selecting a ramp for a public building or public transport, consider all the mobility assistive products available.		Portable ramps		Consider referral for other types of ramp		Other options
		Twin track portable ramp 	Wide platform portable ramp 	Permanent ramp 	Other ramps 	Grab bar/s 
Environment	Entrance to home or personal vehicle	✓	✓	✓	✓	-
	Entrance to public building or transport	✗	✓	✓	✗	-
	Inside home	✗	✓	✓	✓	-
	Inside public building	✗	✓	✓	✓	-
Doorway	Threshold	If person has difficulty crossing threshold discuss other options		✓ Discuss  referral	✓ Discuss  referral	Builder advice on safely lowering/rounding threshold edges
Surfaces	Surfaces uneven and/or not solid in all weather	✗	✗	✓	✗	Discuss ability to change surface
Measurements	Height more than 45 cm	✗	✗	✓	✗	Identify alternative entrances (if available)
	Length more than 350 cm	✗	✗	✓	✗	
Mobility assistive products	Walking aid	✗	✓	✓	✓	Consider grab bars
	Three-wheel wheelchair	✗	✓	✓	✓	-
	Four-wheel wheelchair	✓	✓	✓	✓	-
Lifestyle	Person will be assisted on ramp	✗	✓	✓	✓	-
	Person/assistant cannot carry and position ramp	✗	✗	✓ Discuss  referral	✗	Discuss other options for access
	Person cannot store ramp	✗	✗	✓	✗	